


FlowstarTM

Compact Digital Uroflowmeter

UROLOGY & UROGYNECOLOGY
A practical solution for
everyday uroflowmetry

Precise, cost effective and digital

Flowstar is a compact, portable, precise flowmeter that is designed for practical everyday use to seamlessly integrate with your practice.

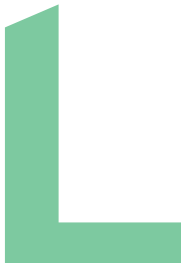

Compact digital uroflowmeter

Trusted and reliable, Flowstar is an optimal uroflowmetry solution for busy clinics where ease of use and precision are priorities.

Flowstar begins recording automatically as soon as the patient starts voiding. Flow/volume scales can be adjusted, and the pre-flow and post-flow times are user selectable.

The integrated printer enables Siroky flow plots to diagnose obstructed voiding.

The lightweight thermal printer unit of the Flowstar is connected by cable and can be up to 10 meters away from the uroflow sensor.

The digital weight sensor can be mounted on a height adjustable uroflow stand, and/or placed under a micturition chair.



©2022 Laborie. All rights reserved.

Flowstar benefits:

- Automatic start-stop feature
- Real-time print-out in portrait or landscape mode
- Reprint of report
- Siroky flow plot
- Adjustable flow and volume scales
- User definable pre- and post-flow time
- Automatic artifact detection software
- High-precision digital weight sensor
- Easy to clean



Laborie

USA

Pease International Tradeport
180 International Drive,
Portsmouth, NH 03801, USA

T +1 802 857 1300
F +1 802 878 1122
W laborie.com



Learn more

103205 - 0